
XTERN AFO – UNIVERSAL LACE-CLIP FITTING

ORTHÈSE XTERN INSTALLATION CLIP-LACET UNIVERSELLE

2022-07-14

-1- SHOE RECOMMENDATIONS


RECOMMENDATIONS SOULIERS

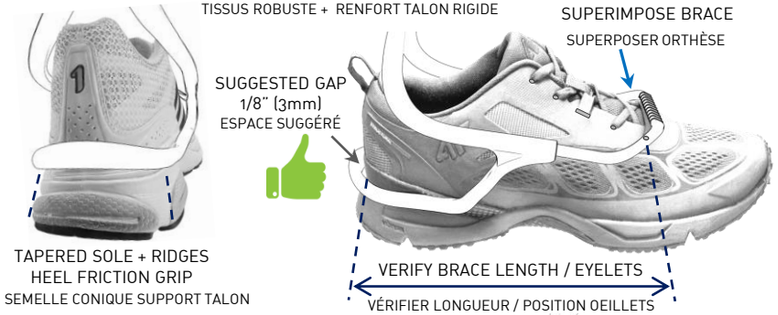

AVOID MINIMALIST SHOES WITH VERY LIGHT FABRIC AND THIN SOLE.
ÉVITER SOULIER MINIMALISTE TISSUS TRES LÉGER ET SEMELLE MINCE.


-2- ANATOMICAL LOCATION FRONT "U"/ LACE CLIP


POSITION ANATOMIQUE "U"AVANT / CLIP-LACET

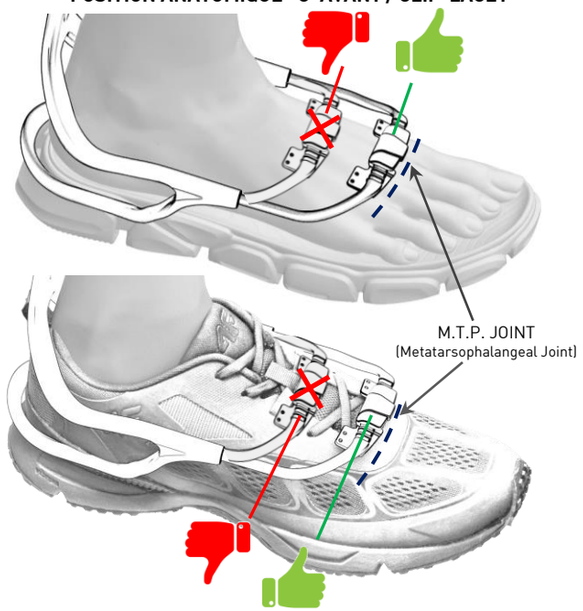




-3-

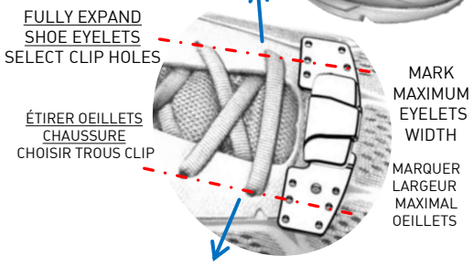


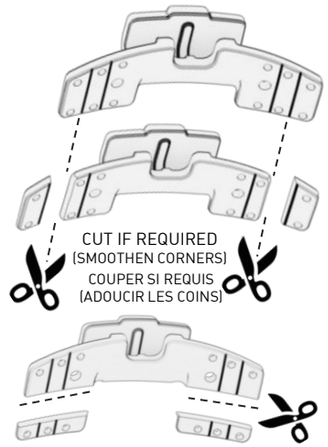


-4-

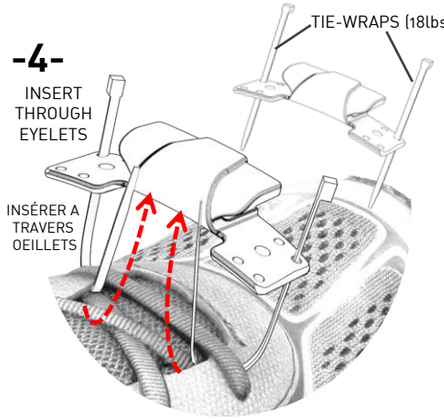


-5-

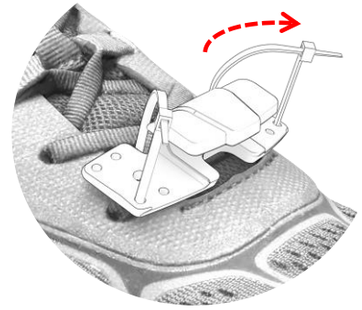


-6-

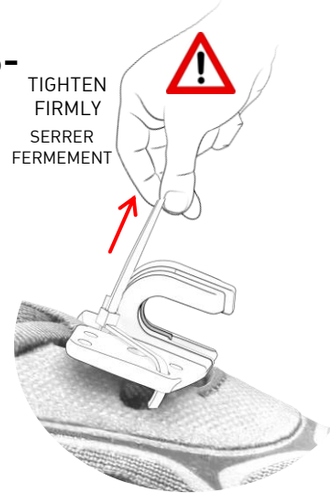

-7-

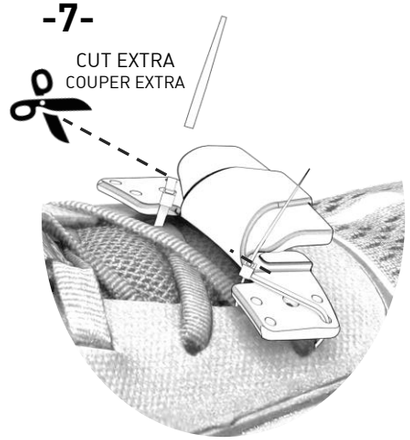

-8-

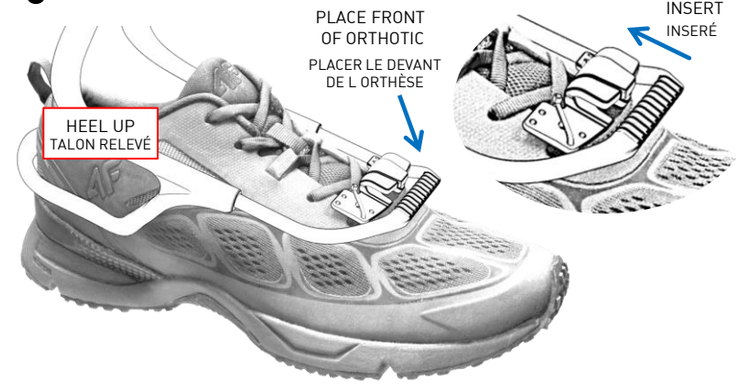

-9-



-10-



1-ENFILER LA CHAUSSURE.

2-ATTACHER LA COURROIE MOLLET/TIBIA.

3-SERRER LES LACETS: LACER SERRÉ AUTOUR DE LA CHEVILLE POUR BIEN SECURISER LE PIED DANS LE SOULIER. GARDER LES LACETS DÉTENDUES PRÈS DU CLIP-LACET.

4-VOUS ETES PRET POUR L ACTION!

1-DON THE SHOE.


2-SECURE THE CALF/SHIN STRAP.

3-TIGHTEN LACES: LACE TIGHT AROUND ANKLE TO SECURE FOOT PROPERLY. KEEP LACES LOOSE NEAR LACE CLIP.


4-YOU'RE READY FOR ACTION!

VELCRO SHOES
SOULIERS À VELCRO

-1-

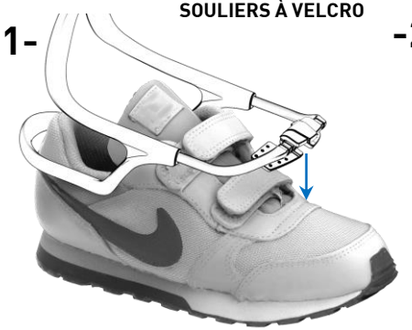


-2-

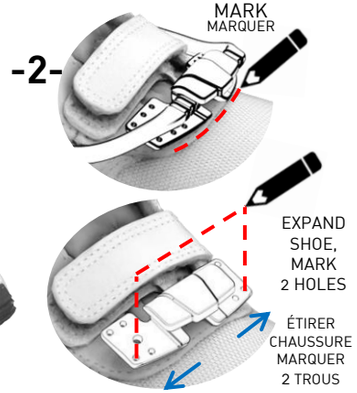


-3-




-4-

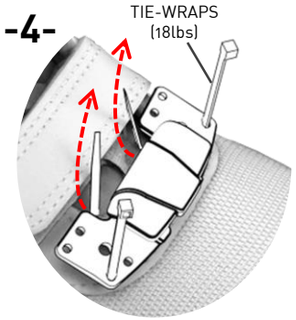


-5-

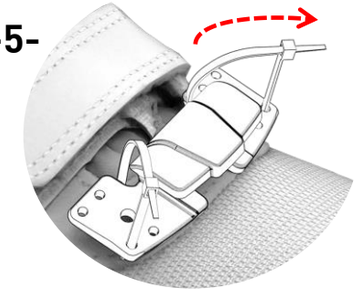

-6- TIGHTEN FIRMLY
SERRER FERMEMENT

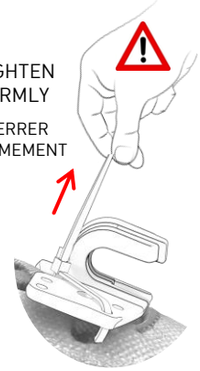


-7-

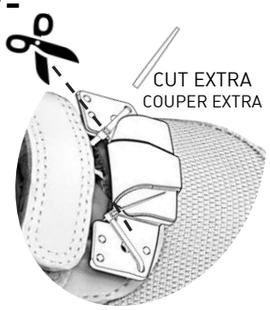




DRESS SHOES
SOULIERS HABILLÉS


-1-

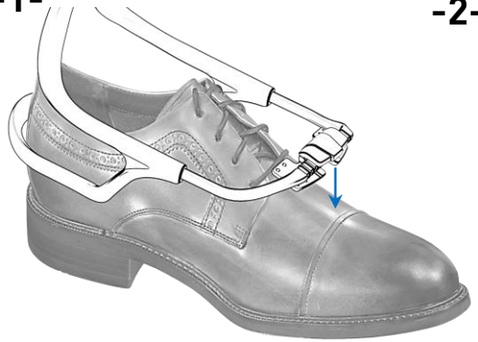


-2-

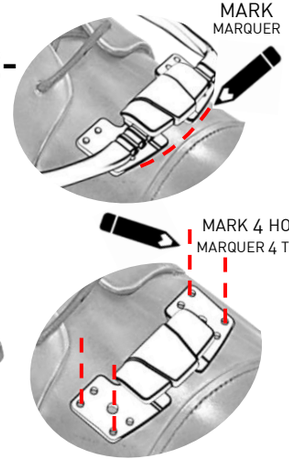


-3-

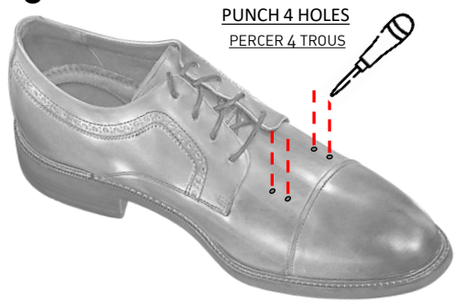


-4-

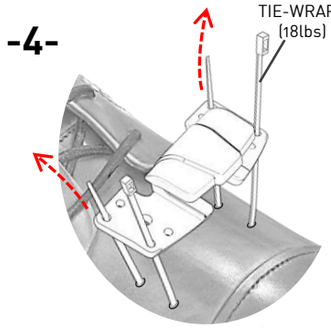


-5-

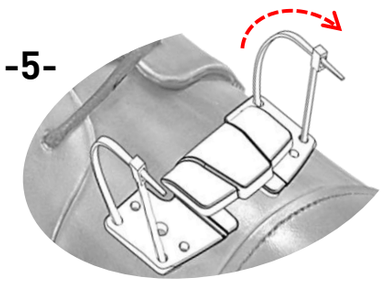

-6- TIGHTEN FIRMLY
SERRER FERMEMENT

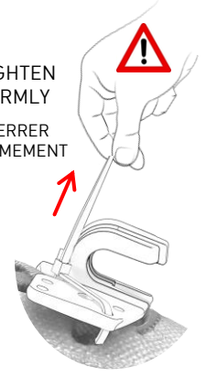


-7-

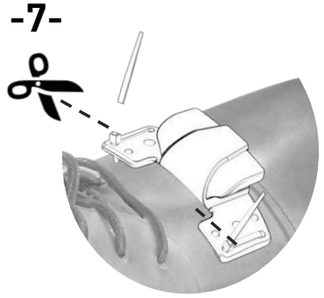



SANDALS
SANDALES


-1-

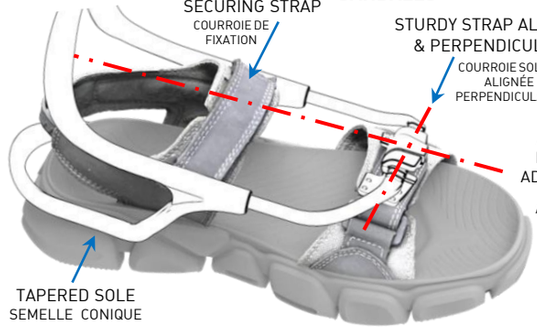

-2-

CHOOSE FIXATION METHOD ACCORDING TO STRAP SITUATION
CHOISIR METHODE FIXATION SELON SITUATION DE LA SANGLE

-3-

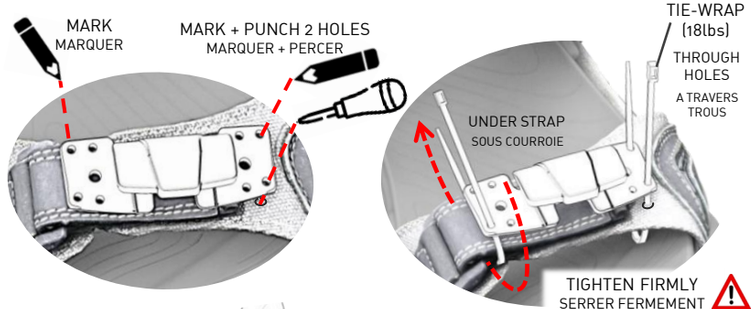



WATER SHOES



IF EYELETS ARE VERY LONG, PUNCH HOLES AT THE BASE OF EYELETS.
SI LES OEILLETS SONT TRÈS LONG, PERCER DES TROUS À LA BASE DES OEILLETS.